
FRED easyport TRAINER

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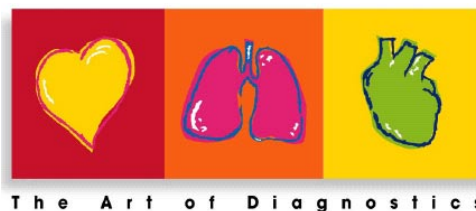Contents

Page

1	General points	19
2	Safety instructions	21
2.1	Responsibilities of the user	21
2.2	Proper use	21
2.3	Organisational measures	21
2.4	Use in accordance with safety rules	21
2.5	Working with other devices	22
2.6	General safety instructions	22
2.7	Other specifications	22
2.7.1	Implicit permission	22
2.7.2	Guarantee conditions	23
2.8	Symbols	23
2.8.1	Symbols used in the manual	23
2.8.2	Symbols used on the device or on the accessories	23
3	Intended use	24
4	Characteristics	24
5	Display and controls	25
5.1	Front views	25
5.1.1	Device	25
5.1.2	Remote control	25
5.2	Rear view device	26
5.3	Electrodes and cable	26
6	Scenario selection procedure	27
6.1	List of available scenarios	27
6.2	Selection of a scenario on the device	27
6.3	Selection of a scenario with the remote control	28
6.3.1	Description of the functions associated with keys on the remote control	28
7	Implementation	29
8	To switch the device off	29
9	Cleaning and maintenance	30
10	Servicing and disposal at the end of the service life	30
10.1	Servicing	30
10.2	Disposal at the end of the service life	30
11	Technical specifications	30
12	Accessories and consumable	31

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1 General points

- The FRED easyport TRAINER device bears the CE mark in accordance with Directive 93/42/EEC relating to medical devices, and fulfils the fundamental requirements of Appendix I of that Directive.
- The device complies with standard EN 55011 class B.
- The safety of the user and trouble free device operation, shall only be guaranteed when original SCHILLER parts are correctly used. Further, only the accessories mentioned in these instructions for use or accessories that have been expressly authorised by SCHILLER may be used. If the device becomes damp as a result of temperature variations, it may only be put back into service once it is completely dry.
- The manufacturer shall not be held liable for the safety, reliability and performance of the device if:
 - Assembly, extensions, adjustments, modifications or repairs are not performed by the manufacturer or by personnel authorised by the manufacturer.
 - The electrical installation of the premises of use does not comply with locally applicable regulations.
 - The device is not used in accordance with the instructions for use.
- The manufacturer may supply, on request, circuit diagrams, lists of components, descriptions, calibration instructions or any other information, for use by the qualified technical personnel of the user for repairing such parts of the device as have been stated to be repairable by the manufacturer. Such supply shall not in any case constitute permission or approval to modify or repair the device.
- This document describes a device type that complies with the safety standards applicable at the time of printing. SCHILLER reserves all rights to the circuits, processes, names, software and devices mentioned in these instructions for use.
- The user guide is an integral part of the device and should always be kept near the device. Close observance of the information given in the user guide is a prerequisite for using the device as intended for and correct operation and ensures operator safety. Therefore, be sure to read the complete user guide.
- The user guide informs the device operator about the intended use, exact function, operation and required preventive maintenance. It is not a substitute for a product training.
- The quality assurance process used by SCHILLER is in compliance with international standards ISO 9001 and ISO 13485.
- No part of these instructions for use may be reproduced without the written permission of SCHILLER.
- This device is not designed for any use other than that specifically described in this manual. Such use may be hazardous.

Manufacturer

SCHILLER Médical SAS
4, rue Louis Pasteur – BP 90050
F-67162 WISSEMBOURG Cedex
Tel: +33 (0)3.88.63.36.00
Fax: +33 (0)3.88.94.12.82
E-mail: info@schiller.fr
Internet site: www.schiller-medical.com

Revision history

Version	Date	Comment
V1.0	November 2006	1st issue

2 Safety instructions

2.1 Responsibilities of the user

- This device may only be used by qualified persons with appropriate training.
- Make sure that the user has read and understood the instructions for use and particularly this "Safety instructions" section.
- Replace any faulty or missing parts immediately.
- The operator shall comply with the applicable regulations in the area of accident prevention and safety instructions.
- The device shall be stored out of the reach of children.
- Dispose of the packaging material and make sure that it remains out of the reach of children.

2.2 Proper use

- Always use the device in accordance with its stated technical specifications.
- Do not use the device in an area with an explosion hazard or in an area with inflammable gases such as anaesthetic gases.

2.3 Organisational measures

- This device may only be used after explanations about its working and safety instructions are given and understood.
- Keep the instructions for use to hand in the place of use of the device. Make sure that they are always complete and legible.
- In addition to the instructions for use, follow the ordinary legal requirements and any other requirements applicable in the area of accident prevention and environment protection.

2.4 Use in accordance with safety rules

- Read and follow the instructions for use, particularly the safety instructions.
- Immediately report to the person responsible any change liable to affect safety, including changes relating to the behaviour of the device in operation.
- Ensure of the correct operation of the device before each use
- Only connect original SCHILLER accessories to the device.

- Before starting the device, check that there is no damage to its housings.
- The device may only be exposed to wide temperature variations for a very short time. Very high temperature variations may lead to the formation of condensation water on the device. However, if condensation water does form, dry the device and all its connections.
- Install the device so that it cannot fall on the user or the floor.
- The device must not be connected nor used on a patient.

2.5 Working with other devices

- Use only the accessories supplied and recommended by SCHILLER. The use of accessories manufactured by third parties may lead to injuries, provide inaccurate information and/or lead to device damage.
- Very high leakage currents (total leakage currents) can be a hazard for the user.
- The working of the device may be influenced by the presence of magnetic or electrical fields from X-ray machines, tomography machines, radio transmitters and mobile telephones, among other devices. Avoid the use of such devices and keep a sufficient distance from them.

2.6 General safety instructions

- Risk of injury for the user if the device is used with faulty housings or cables! That is why the following instruction must be followed strictly:
Immediately replace any faulty device, cable or connector.

2.7 Other specifications

2.7.1 Implicit permission

- The possession or acquisition of this device does not amount to express or implicit permission to use this device with parts which, alone or in combination with this device, may be covered by one or more patents.

2.7.2 Guarantee conditions

Your **FRED easyport TRAINER** device from SCHILLER comes with a guarantee (from the date of purchase) covering all workmanship or material defects. The guarantee does not however cover damage due to negligence or improper use. The guarantee includes the replacement of faulty parts. It excludes any liability for consequential damage. The guarantee shall not apply if the device is repaired by unauthorised or unqualified personnel.

In the event of a fault, the device in question shall be returned to the local SCHILLER representative or directly to the manufacturer. The manufacturer warrants the safety, reliability and proper operating of the device only if

- any assemblies, additions, new settings, modifications or repairs are carried out by persons approved by the manufacturer and
- the device and its accessories are used in accordance with the manufacturer's specifications.

Important

SCHILLER offers no other guarantee relating to the merchantable quality or the fitness of the product or its components for a particular purpose.

2.8 Symbols

2.8.1 Symbols used in the manual

The device meets the essential requirements of directive 93/42.

2.8.2 Symbols used on the device or on the accessories

The device meets the essential requirements of directive 93/42.

Do not dispose with household waste. Dispose separately.

Caution : refer to accompanying documentation.

3 Intended use

- **FRED easyport TRAINER** is a simple cost-effective simulator of the FRED easyport that is designed to be used exclusively for training.
- **In no event FRED easyport TRAINER may be used as a conventional defibrillator. It should not be connected to a patient.** It does not deliver any energy.
- **FRED easyport TRAINER** should only be used on a mannequin with specific training electrodes.
- **FRED easyport TRAINER** has been designed for the purpose of training users of SCHILLER's FRED easyport automated defibrillator. The device provides several simulations or scenarios (from 0 to 8) that aim to help to familiarise with the FRED easyport and shows that the users have the knowledge required to use it in case of emergency. The device also allows the trainer to work using an infrared remote control to select a scenario or to intervene in real time on the course of an intervention.

4 Characteristics


- **FRED easyport TRAINER** allows the **simulated** delivery of electric shocks. There is no high voltage in the device, and training is risk free.
- **FRED easyport TRAINER** may be used on all types of mannequins as it does not deliver an electric shock. It may only be used with the training electrodes sold by Schiller.
- **FRED easyport TRAINER** may be used with an infrared remote control. The remote control allows the trainer to control the training actions.
- **FRED easyport TRAINER** has 9 pre-recorded training scenarios which simulate realistic cardiac arrest configurations.

The basic equipment of the FRED easy TRAINER consists of :

- The device **FRED easyport TRAINER** shipped with one 9 Volts alkaline battery
- The remote control RC-7 shipped with two 1,5 Volts alkaline batteries
- Single use training electrodes
- The reusable connecting cable for training. The connector of the cable is fitted with a magnet allowing the device to recognise the presence of the electrodes without any electrical connection.

5 Display and controls

5.1 Front views

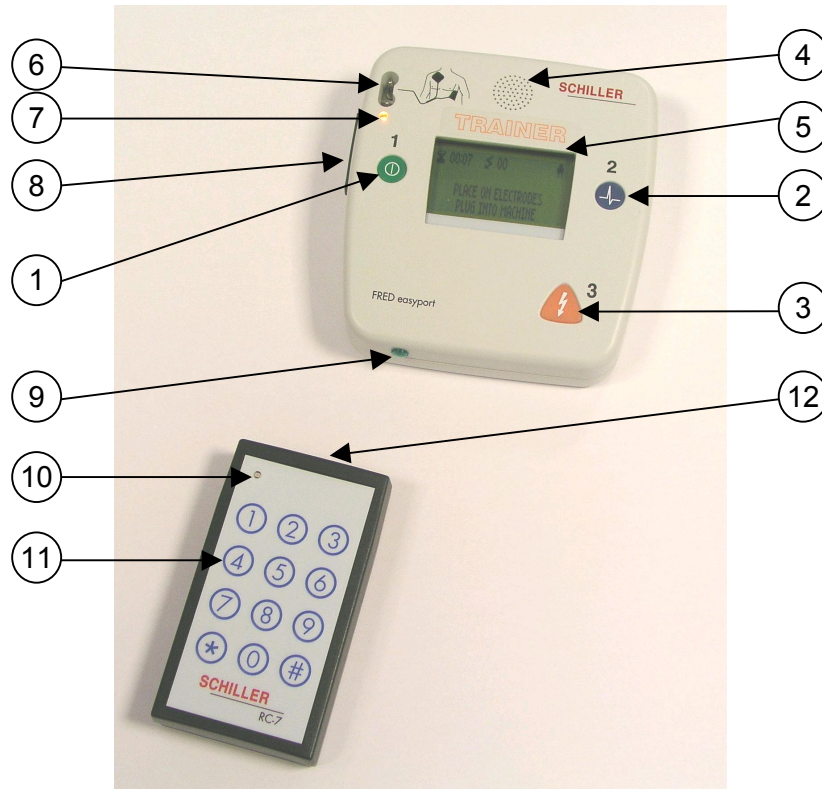
5.1.1 Device

- 1 Green On/Off button
- 2 Blue ANALYSE button
- 3 Orange SHOCK button for shock simulation
- 4 Loud speaker
- 5 LCD display : displaying of elapsed time, number of shocks simulated, low battery and text messages
- 6 Connector for connecting cable of the training electrodes
- 7 Diode for electrode fault simulation
- 8 Connector for voice messages downloading in case of language changing or message modification
- 9 Infrared receiver of data emitted by the remote control

5.1.2 Remote control

- 10 Control diode (lights when pressing on a key)
- 11 Keyboard for scenario selection
- 12 Infrared emission area


5.2 Rear view device

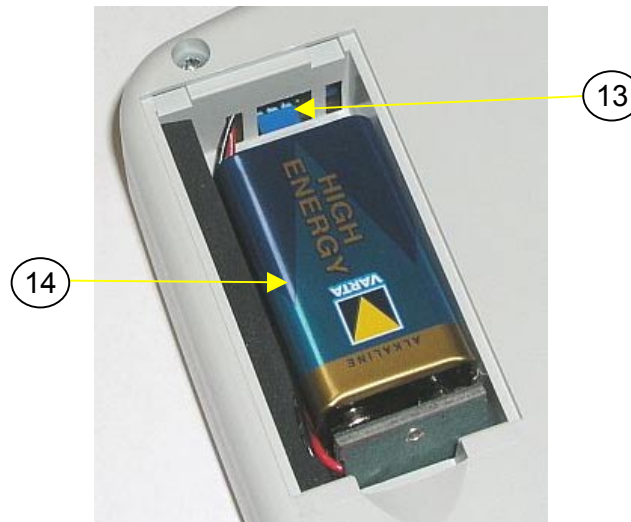13 Connector for downloading a new software

14 9 Volts alkaline battery

5.3 Electrodes and cable

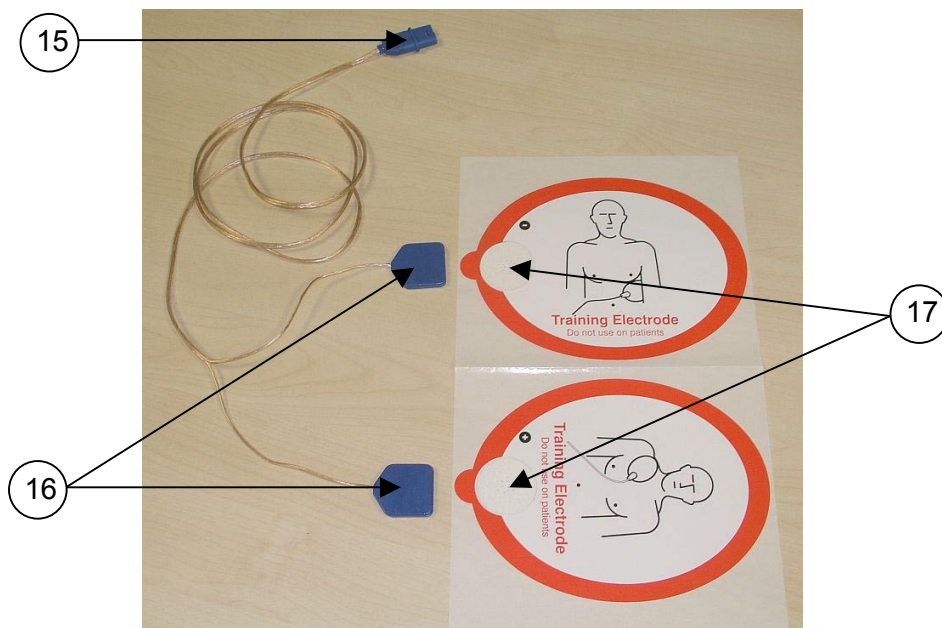15 Connector of the connecting cable of the electrodes fitted with a magnet allowing the device to simulate a good contact of the presence on the mannequin without electrical contact

16 Ends of the cable to be stuck on the electrodes

17 Self adhesive discs of the training electrodes

6 Scenario selection procedure

6.1 List of available scenarios

Scenario 0	SNR
Scenario 1	SNR / SR
Scenario 2	SNR / SR / SNR
Scenario 3	SR
Scenario 4	SR / SNR
Scenario 5	SR / SNR / SR
Scenario 6	Elect fault / SR
Scenario 7	SR/ Elect fault / SR
Scenario 8	SR / SR / SNR
Scenario 9	Infrared control function

Notes :

SNR	Shock not recommended
SR	Shock recommended
Elect fault	Electrode fault

The last step of each scenario is indefinitely repeated.

6.2 Selection of a scenario on the device

- Keep the orange SHOCK button (3) pressed while putting the device On with the green On/Off button (1) to access in the configuration menu.
- The selection in a menu is made by the blue ANALYSE button (2); the passage to the next step is made by the orange SHOCK button (3).
 - Line 1: functioning mode of the device
 - Line 2: choice between automatic and manual analysis
 - Line 3: selection of one of the scenarios (0 to 8) or selection of the remote control (scenario 9).
 - Line 4: saving of the selected configuration with the blue ANALYSE button (2): the device goes Off.

Note: pressing longer than 3 seconds on the green On/Off button (1) puts the device Off without saving the settings.

The device can then be switched on normally and it will operate according to the selected scenario.

6.3 Selection of a scenario with the remote control

The trainer has the possibility to select a scenario with the remote control when the device On. Press on # key, slacken it, then press on the key corresponding to the selected scenario; the indicator (10) of the remote control lights when a key is pressed.

Notes:

- When the choice of the scenario is made during the initialization of the device (displaying of the message "FRED easyport is testing"), this one is immediately taken into account.
- When the choice of the scenario is made after the initialization of the device, this one is taken into account with the next starting of the device.

In the case of the selection of scenario number 9, the device waits for the remote control's orders to change it's status. If the device don't receive orders from the remote control, the shock will always be recommended like in scenario 3.

6.3.1 Description of the functions associated with keys on the remote control

Key 0	Immediate move to the next step
Key 1	Restart at the analysis step; reset of the counters of time and number of delivered shocks
Key 2	Shock recommended, during the analysis in progress or during the next analysis
Key 3	Shock not recommended, during the analysis in progress or during the next analysis
Key 4	The rhythm has changed. This message appears after a recommendation of shock
Key 5	Press on the electrodes (bad contact). This message appears during analysis
Key 6 *	Movements detected. Stop movements. This message appears during analysis
Key 7	Extinction of the device
Key 8	Resume CPR. This message is immediately taken in account.
Key 9	Pause : the current scenario is interrupted; the display is frozen. Press again on key 9 to continue the training.

* This function is not available on the defibrillator FRED easyport

7 Implementation

- Remove the self-adhesive disks (17) from the training electrodes.
- Stick the blue ends (16) of the connecting cable on the electrodes.
- Begin the training

Note:

Use of the remote control:

- The transmitter of the remote control (12) must be directed in direction of the receiver (9) on FRED easyport TRAINER.
- No object or nobody must be between the remote control and the receiver.
- The angle of transmission between the infra-red beam and the transmitter or the receiver should not exceed 45°.

8 To switch the device off

Press longer than 3 sec on the green On/Off button (1) to switch the device off.

9 Cleaning and maintenance

The outer surfaces of the device may be cleaned with cotton wool moistened with alcohol.

CAUTION ! Risk of damage to the device. Never use phenol based cleaning agents or agents containing peroxide derivatives to disinfect the surfaces of the device housing.

During cleaning, the batteries must be removed from the device. No liquid may be allowed to enter the device. If liquid does penetrate into the device, clean the device and have it verified completely.

10 Servicing and disposal at the end of the service life

10.1 Servicing

Check the device performance once a year.

10.2 Disposal at the end of the service life

At the end of its service life, the device and the accessories must be disposed of in compliance with the local regulations.

The device must be disposed in a municipally approved collection point or recycling centre when it is no longer required.

Improper disposal harms the environment and human health due to the presence of dangerous substances in electrical and electronic equipment.

11 Technical specifications

- Operating mode : training in semi-automatic defibrillation
- Detection of electrode presence by means of a magnet placed in the connector of the connecting cable for the training electrodes
- **Power supply of FRED easyport TRAINER**
 - 1 alkaline battery type 6LR61 9 Volts
 - Displaying of a symbol when battery low
- **Infrared remote control**
 - Power supply : 2 alkaline batteries type LR6 1.5 Volts
 - Range : 4 m
 - Maximum angle of deviation compared to the transmitter or the receiver : 45°

Physical and environment requirements

- Transport / storage:
Temperature : – 20 to + 50 °C
Relative humidity rate : 30 to 95 %, without condensation
Atmospheric pressure : 500...1060 hPa
- Operation:
Temperature : 0 to + 40 °C
Relative humidity rate : 30 to 95 %, without condensation
Atmospheric pressure : 700 to 1060 hPa
- Protection case : IPX 4

Dimensions and weight

Width	125 mm
Depth	130 mm
Height	37 mm
Weight	< 300 g with battery

12 Accessories and consumable

Part Number	Designation
1-67-9900	FRED easyport TRAINER
1-67-5400	Remote control
0-21-0019	Reusable cable for training electrodes ref. 0-21-0005. Length: 1.40 m
0-21-0005	Adhesive training electrodes (10 sets of 2 electrodes)
35492	9 Volt alkaline battery for FRED easyport TRAINER
53108	1.5 Volt alkaline battery for remote control (1 box of 4 batteries)

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